

Ipas[®] M·V·A

MANUAL VACUUM ASPIRATORS



product catalog

IpasMVA.com

compassionate healthcare

WELLNESS CREATED THROUGH EDUCATION AND EMPATHY

The Ipas MVA was thoughtfully designed to be a portable, hand-held device for uterine evacuation. Unlike the complex and loud electric vacuum aspirators, our MVA is easily used by a provider in a regular exam room or ER setting. The benefit to providers and facilities is a reduction of cost and time—which means you can provide better care for more patients.



15-minute.....
UTERINE EVACUATION

97-99% SUCCESS RATE

10-15 MINUTE PROCEDURE

NO SEDATION REQUIRED

MULTI-SETTING VERSATILITY

LOW RISK OF COMPLICATIONS

WORLDWIDE USAGE

OVER 100 COUNTRIES • OVER 50 YEARS

The Ipas Manual Vacuum Aspirator revolutionized access to uterine evacuation by reducing both time and expense. Easy to use and simple to learn, it gained rapid use from providers eager

to bypass operating theaters due to its suitability for outpatient settings. It continues to grow in worldwide usage due to its efficacy and safety.



Ipas MVA PLUS® - Reusable up to 25 times

“

MVA is a safe, reliable and effective method of evacuation of the uterus in early pregnancy, less than 12 weeks... We have replaced conventional D&E, D&C and suction curettage... with Manual Vacuum Aspiration.

– Prof. Ghazala Mahmud, Head of Department of Obstetrics and Gynaecology, Fazia Medical College

patient-centric care

EVERY WOMAN'S INTELLIGENT AND NURTURING SOLUTION

Every minute matters. The Ipas® MVA has been proven to provide safe and effective uterine evacuation in under 15 minutes. Because the device can be used in an outpatient or exam-room setting, it also offers peace of mind to your patients. They can return to their lives, work, or families the same day, and recover in the comfort of home. Choosing our Ipas MVA is a reliable way to offer both health and hope. Simultaneously.

A BENEFIT FOR PATIENTS & FAMILIES



RELIABLE & SAFE

SHORT RECOVERY

QUIET ~ NO MACHINES

FEW SIDE EFFECTS*



*Minimal pain/cramping



PROVIDING PATIENTS
recovery options
THAT MATTER TO THEM



our benefits & usage

THE MVA — A GENTLE APPROACH FOR ALL PROVIDERS

The Ipas MVA—a hand-held vacuum that doesn't require electricity—uses suction to remove uterine contents through the cervix using a cannula. Easy to learn and use by nurse practitioners, midwives, and other primary healthcare providers, the MVA procedure is safe, 98% effective, affordable, and uses minimal resources.




Note: The diagram above showcases our Double Valve. All these attributes also apply to the Ipas MVA PLUS®.

PAIN MANAGEMENT PHARMACOLOGICAL OPTIONS

There are several options for anesthesia or sedation for MVA. Discussed with the patient, the chosen method should also be relevant results of a clinical assessment. Methods may include:

- Nonsteroidal anti-inflammatory drugs
- Local anesthetic injection
- Monitored anesthetic care (MAC)

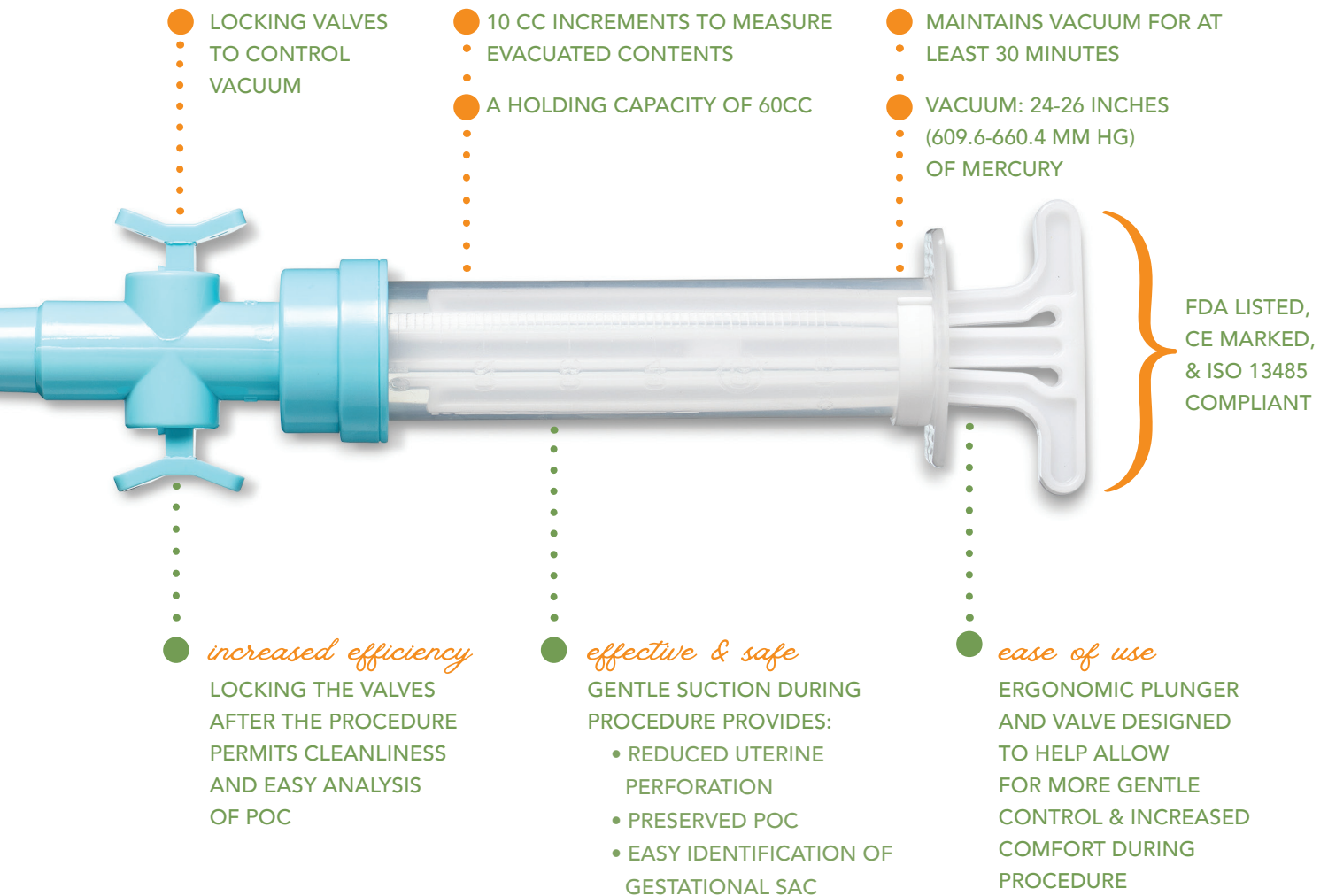

EDUCATION & ENVIRONMENT

Pain relief techniques can be combined with pharmacological interventions. Education, distraction techniques, controlled breathing, and

listening contribute to well-being. The multi-setting versatility means the MVA can be used in a less clinical and more nurturing exam room/setting.

Ipas[®] MVA KITS


Item Code: UEK

MISCARRIAGE MANAGEMENT KITS

An MVA ideal for emergency departments, ambulatory surgery centers, and office-based settings. Our single-use kits are pre-packaged, simple, and portable.

EACH KIT INCLUDES:

- 1 Ipas Double Valve Aspirator
- 4 Ipas EasyGrip[®] Cannulae
(6, 7, 8, 9mm)
- 3 Denniston Dilators
(5/6, 7/8, 9/10 mm)

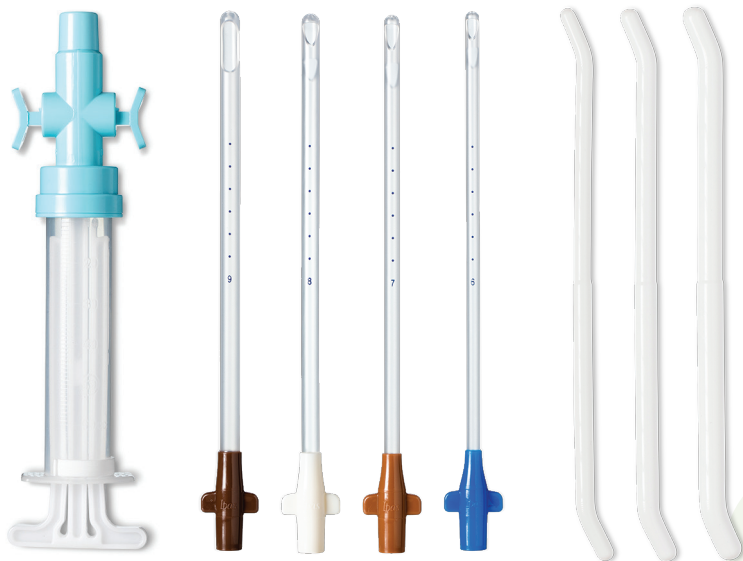

“

I use MVA exclusively for all office-based aspiration procedures... The reusable Plus or disposable kit makes it convenient for the ER, clinic, or office. It's quick, easy to use, and my patients appreciate having another option.

– Mitchell Crenin, MD
Director of Family Planning, University
of California, Davis Health System





ENDORSEMENTS FROM THE INDUSTRY

Vacuum aspiration (MVA and EVA) is recommended over D&C by the World Health Organization and the International Federation of Gynecology and Obstetrics. Most major professional obstetric and

gynecological groups, like the American College of Obstetricians and Gynecologists (ACOG) and the Royal College of Obstetricians & Gynaecologists (RCOG), have issued similar position statements.

Ipas[®] ASPIRATORS

EASY TO LEARN & SIMPLE TO USE

Regardless of which aspirator you choose providers can quickly learn to use the device with minimal training.

DOUBLE VALVE MVA

EFFECTIVE, AFFORDABLE & VERSATILE

For customers who prefer the simplicity of single-use devices. A low-cost and flexible option.

- Used In all 50 states for more than 25 years
- Compatible with cannulae sizes 4-10mm and 12mm
- Individually wrapped, sterile, and ready to use



Item Code: DVSSU

**THE ONLY
DOUBLE VALVE MVA
MANUFACTURED IN
A USFDA INSPECTED
AND REGISTERED
FACILITY**

MVA PLUS[®]

The Ipas MVA Plus is the most advanced MVA in the world:

- Designed for reuse
- Can be disassembled for easy cleaning
- Compatible with cannulae sizes 4-10mm and 12mm
- The only MVA in the world that tolerates heated methods (i.e., autoclave or boiling water) for processing.



Item Code: PLUS

Both aspirators are FDA listed, CE marked, and ISO 13485 Compliant



“

We use the MVA in our hospital routinely for outpatient procedures. All postgraduate students in obstetrics and gynecology should learn this technique. All practicing gynecologists should adopt this method of uterine evacuation.

*– Shershah Syed, FRCOG,
Pakistani physician and surgeon*

OUR IPAS ASPIRATORS HAVE:

- A holding capacity of 60CC
- 10 CC increments to measure evacuated contents
- Locking valves to control vacuum
- Ability to maintain vacuum for at least 30 minutes
- Vacuum: 24-26 inches (609.6-660.4 MM HG) of mercury
- Shelf life is 5 years from month of manufacture

Ipas[®] CANNULAE

SEMI-RIGID & FLEXIBLE OPTIONS


Designed specifically for use and integration with the Ipas MVA, our cannulae are available in multiple sizes.

IPAS EASYGRIP[®]

Eight color-coded cannulae. The tactile responses of a rigid curette with the gentle probe of a flexible cannula. Used with the Ipas MVA Plus[®] or Ipas Double Valve Aspirator.

- Semi-rigid
- No adapters required

Item Code: SR#SU

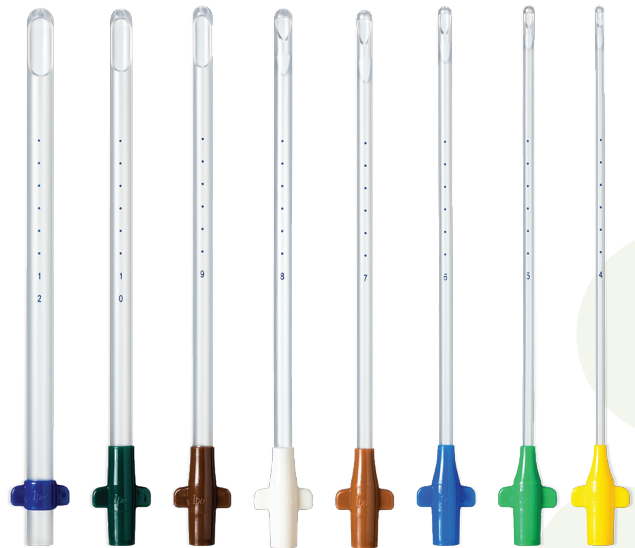

FLEXIBLE KARMAN

Softer, more flexible cannulae may be a better option for retroverted or retroflexed uterine position. Sizes 4-6mm fit the Ipas Single Valve Aspirator.

- Adapter required if used with Ipas MVA Plus[®] or Ipas Double Valve

Item Code: FU#SU



ADAPTER SET

For use with the Flexible Karman



Item Code: ADPSET

To order any cannula separately, replace item code # with diameter size.



ALL OUR CANNULAE ARE:

- Sterile & individually packaged
- Disposable
- Color coded by size
- Ready to use
- Available in 8 diameter sizes (4, 5, 6, 7, 8, 9, 10, and 12mm)

Ipas[®] DEVICES

PRODUCTS & OTHER DEVICES

DENNISTON DILATORS

Reusable, steam autoclavable at 250°F, and processable by cold methods — Cidex[®] or Sporax[®] II.

EACH SET INCLUDES:

5 double-ended, tapered dilators
(5/6, 7/8 9/10, 11/12, 13/14mm)



Item Code: DD

PLUS ACCESSORIES

Replacement parts & silicone for use with the PLUS aspirator

Item Code: AKPL

Code: SILICONE



SINGLE VALVE

Designed to work with the 3mm cannula, this single valve aspirator is ideal for use in Endometrial Biopsy procedures.

- Smaller device
- Ideal for Menstrual Regulation
- For procedures up to 6 gestational weeks
- Available sterile for contamination-free sampling with endometrial biopsy 3mm cannula
- Compatible with 4-6mm Flexible Karman cannulae



Item Code: SVSSU

3MM CANNULA

Designed for use with the Ipas Single Use Aspirator.



Item Code: EB3

LET'S PURSUE

health & comfort

SIMULTANEOUSLY...

BECAUSE A FEW

quiet minutes

CAN GRANT WELLNESS

**TO ORDER
THE IPAS MVA:**

(800) 850-1657
or
customerservice@hpsrx.com
HPSRx.COM



We also provide:

Dilapan-S
Med/Surg
Pharmaceuticals
Birth Control
24/7 Online Ordering
US Customer Service
Shipping to 50 States

For more details:
IpasMVA.com